


INSTRUCTIONS FOR USE OF

KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

August 2025, Version 7.0

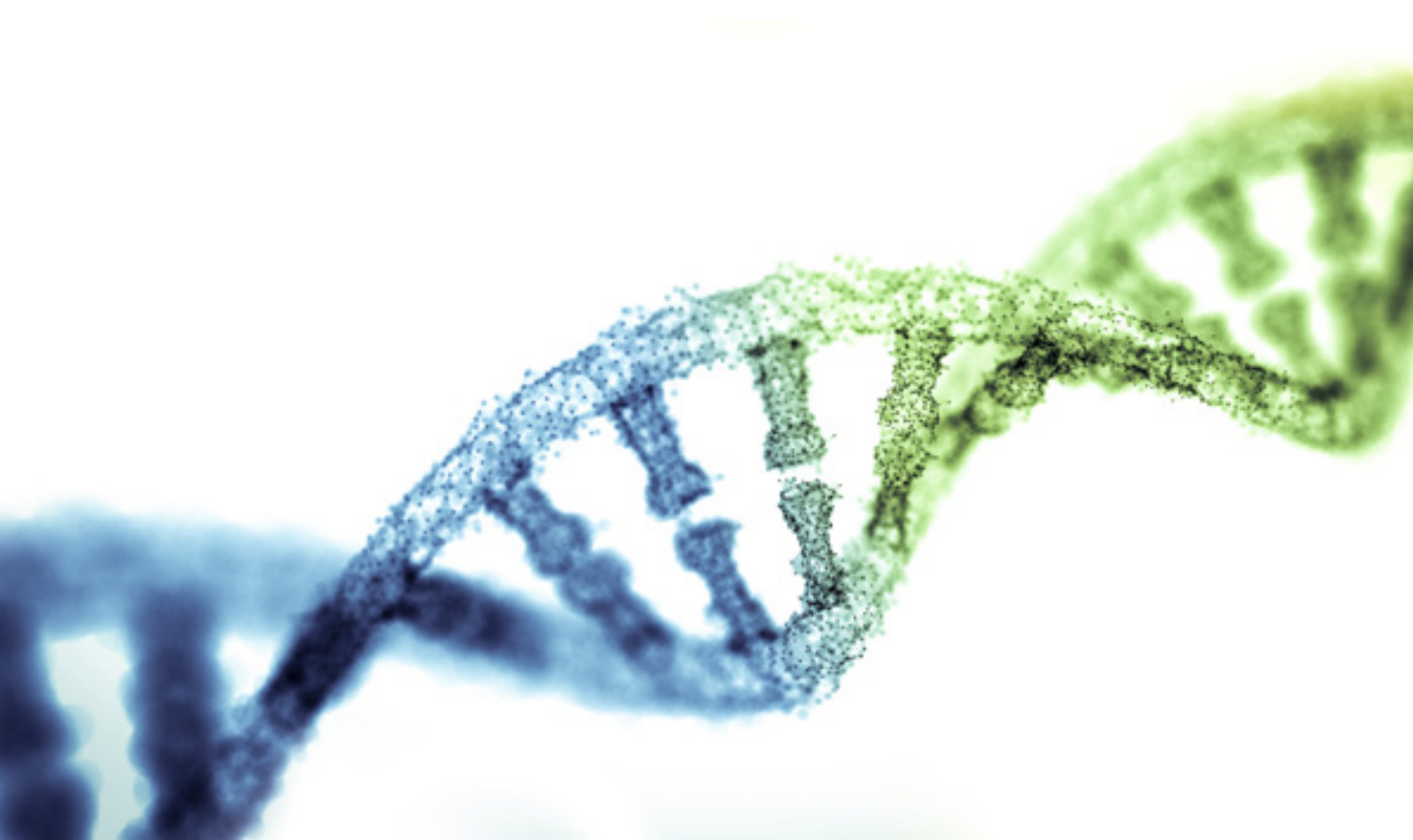

Contents

The KAPA RNA HyperPrep Kits with RiboErase (HMR) Globin contain:

Kit Material Number	Description	Volume
KK8562 08308314702 24 Libraries consists of the following sub kits: KK8542 08105936001 and KK8546 08105863001 and KK8481 07962266001 and KK8478 08308276001	Hybridization Buffer	110 µL
	Hybridization Oligos (HMR)	110 µL
	Globin Hybridization Oligos (HMR)	24 µL
	Depletion Buffer	80 µL
	RNase H	55 µL
	DNase Buffer	60 µL
	DNase	55 µL
	Fragment, Prime and Elute Buffer (2X)	300 µL
	1st Strand Synthesis Buffer	300 µL
	KAPA Script	25 µL
	2nd Strand Marking Buffer	780 µL
	2nd Strand Synthesis & A-Tailing Enzyme Mix	50 µL
	Ligation Buffer	1 mL
	DNA Ligase	280 µL
	PEG/NaCl Solution	1 mL
KAPA Pure Beads (2 x 3.2 mL bottles)	6.4 mL	
Library Amplification Primer Mix (10X)	138 µL	
KAPA HiFi HotStart ReadyMix (2X)	690 µL	

Kit Material Number	Description	Volume
KK8563 08308241702 96 Libraries consists of the following sub kits: KK8544 08105952001 and KK8001 07983280001 and KK8482 07962274001 and KK8479 08308284001	Hybridization Buffer	480 µL
	Hybridization Oligos (HMR)	480 µL
	Globin Hybridization Oligos (HMR)	96 µL
	Depletion Buffer	360 µL
	RNase H	240 µL
	DNase Buffer	264 µL
	DNase	240 µL
	Fragment, Prime and Elute Buffer (2X)	1.4 mL
	1st Strand Synthesis Buffer	1.4 mL
	KAPA Script	130 µL
	2nd Strand Marking Buffer	3.8 mL
	2nd Strand Synthesis & A-Tailing Enzyme Mix	250 µL
	Ligation Buffer	5 mL
	DNA Ligase	1.26 mL
PEG/NaCl Solution	5 mL	
KAPA Pure Beads	30 mL	
Library Amplification Primer Mix (10X)	600 µL	
KAPA HiFi HotStart ReadyMix (2X)	3 mL	

Note: Accessory kits (KAPA Universal Adapter & KAPA UDI Primer Mixes and KAPA Unique Dual-Indexed Adapters) are sold separately.

Shipping, Storage and Stability

- KAPA Pure Beads are shipped on dry ice or ice packs, depending on the destination country. Upon receipt, store KAPA Pure Beads at +2°C to +8°C.
- Enzymes and buffers for rRNA & globin depletion, cDNA synthesis and library preparation are shipped on dry ice or ice packs, depending on the destination country. These components are temperature sensitive, and appropriate care should be taken during storage.
- Upon receipt, store these reagents at -15°C to -25°C in a constant-temperature freezer.
- The 1st Strand Synthesis Buffer, PEG/NaCl Solution, and KAPA Pure Beads are light sensitive and should be protected from light during storage.
- When stored under these conditions and handled correctly, the kit components will retain full activity until the expiry date indicated on the kit label.

Applications

The KAPA RNA HyperPrep Kit with RiboErase (HMR or Human/Mouse/Rat) Globin contains all of the buffers and enzymes required for depletion of rRNA and globin mRNA transcripts and the construction of stranded RNA-Seq libraries from 25 ng – 1 µg of purified, blood-derived total RNA. The KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin is designed for both manual and automated NGS library construction. The kit depletes both cytoplasmic (5S, 5.8S, 18S, and 28S), and mitochondrial (12S and 16S) rRNA species as well as globin mRNA transcripts. The protocol is applicable to a wide range of RNA-Seq applications, including:

- gene expression analysis of high- and low-quality RNA samples (e.g., extracted from FFPE tissue)
- single nucleotide variation (SNV) discovery
- splice junction and gene fusion identification
- characterization of both polyadenylated and nonpolyadenylated RNAs, including noncoding and immature RNAs.



This kit is not compatible with small RNAs <100 bp in length.

Warnings and Precautions

- Wear the appropriate personal protective equipment, such as gloves and safety glasses, to avoid direct contact while handling the reagents. Wash hands thoroughly after handling samples and reagents.
- Use good laboratory practices to avoid contamination when working with the reagents.
- Do not eat, drink, or smoke in the laboratory area.
- Do not pipette by mouth.
- In the event of a spill, clean up the solution with absorbent pads, allow pads to dry, and dispose of pads.
- Observe all national, regional, and local regulations for waste disposal and management.
- Safety Data Sheets (SDS) are available [online](#) or upon request from the local Roche office.

Ordering Information

For a complete overview of Roche Sequencing products, including KAPA RNA HyperPrep Kits with RiboErase (HMR) Globin, go to sequencing.roche.com products.

Trademarks

KAPA and KAPA HYPERPLEX are trademarks of Roche.

All other product names and trademarks are the property of their respective owners.

Contact and Support

If you have questions, contact your local Roche Technical Support. Go to sequencing.roche.com/support for contact information.

Copyright

© 2015 – 2025 Roche Sequencing Solutions, Inc. All Rights Reserved.

Roche Sequencing Solutions, Inc.

4300 Hacienda Drive

Pleasanton, CA 94588 USA

Editions

Version 1.17, June 2017; Version 2.17, September 2017; Version 2.17, October 2018 (republished); 3.20, December 2020; Version 4.21, April 2021; Version 5.23, June 2023; Version 6.0, February 2024; Version 7.0, August 2025

Changes to Previous Versions

Minor updates to align with latest product references and specifications.



Restrictions and Liabilities

This document is provided “as is” and Roche Sequencing Solutions, Inc. (Roche) assumes no responsibility for any typographical, technical, or other inaccuracies in this document. Roche reserves the right to periodically change information that is contained in this document; however, Roche makes no commitment to provide any such changes, updates, enhancements, or other additions to this document to you in a timely manner or at all.

OTHER THAN THE LIMITED WARRANTY CONTAINED IN THIS USER GUIDE, ROCHE MAKES NO REPRESENTATIONS, WARRANTIES, CONDITIONS OR COVENANTS, EITHER EXPRESS OR IMPLIED (INCLUDING WITHOUT LIMITATION, ANY EXPRESS OR IMPLIED WARRANTIES OR CONDITIONS OF FITNESS FOR A PARTICULAR PURPOSE, NON-INFRINGEMENT, MERCHANTABILITY, DURABILITY, TITLE, OR RELATED TO THE PERFORMANCE OR NON-PERFORMANCE OF ANY PRODUCT REFERENCED HEREIN OR PERFORMANCE OF ANY SERVICES REFERENCED HEREIN).

This document might contain references to third party sources of information, hardware or software, products, or services and/or third party web sites (collectively the “Third-Party Information”). Roche does not control, and is not responsible for any Third-Party Information, including, without limitation the content, accuracy, copyright compliance, compatibility, performance, trustworthiness, legality, decency, links, or any other aspect of Third-Party Information. The inclusion of Third-Party Information in this document does not imply endorsement by Roche of the Third-Party Information or the third party in any way.

Roche does not in any way guarantee or represent that you will obtain satisfactory results from using Roche products as described herein. The only warranties provided to you are included in the Limited Warranty enclosed with this guide. You assume all risk in connection with your use of Roche products.

Roche is not responsible nor will be liable in any way for your use of any software or equipment that is not supplied by Roche in connection with your use of Roche products

Conditions of Use

You are responsible for understanding and performing the protocols described within. Roche does not guarantee any results you may achieve. These protocols are provided as Roche’s recommendations based on its use and experience with Roche products.

Use Restrictions

For patent license limitations for individual products, refer to: www.technical-support.roche.com.

01

Instructions for Use 1

Contents	2
Shipping, Storage and Stability	3
Applications	4
Warnings and Precautions	4
Ordering Information	4
Trademarks	4
Contact and Support	4
Copyright	4
Editions	4
Changes to Previous Versions	4
Restrictions and Liabilities	5
Conditions of Use	5
Use Restrictions	5

02

Preface 9

Regulatory Disclaimer	9
Contact Information	9
Conventions Used in This Manual	10

03

Chapter 1. Before You Begin 11

Figure 1 - KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin Quick Guide	12
Protocol Information & Safety	13
References	13
Terminology	14
Prepare the Following Equipment and Reagents	14
Required Equipment, Labware & Consumables	15

04

Chapter 2. Store and Prepare the Reagents 17

Step 1. Store the Reagent Kits	18
Step 2. Prepare KAPA Adapters and KAPA UDI Primer Mixes	18

05

Chapter 3. Prepare the Sample Library 21

Sample Requirements	23
Step 1. Oligo Hybridization and rRNA /globin mRNA Depletion	24
Step 2. rRNA /Globin mRNA Depletion Cleanup	25
Step 3. DNase Digestion	26
Step 4. DNase Digestion Cleanup	27
Step 5. RNA Elution, Fragmentation and Priming	28
Step 6. 1st Strand Synthesis	29
Step 7. 2nd Strand Synthesis & A-tailing	30
Step 8. Adapter Ligation	31
Step 9. Purify the Sample Library using KAPA Pure Beads	33
Step 10. Purify the Sample Library using PEG/NaCl	34

06

Chapter 4. Amplify the Sample Library 35

Step 1. Prepare the Library Amplification Reaction	37
Step 2. Perform the Library Amplification	38
Step 3. Purify the Amplified Sample Library using KAPA Pure Beads	39

07

Chapter 5. Quality Control 41

Quantification	42
Sizing	42

08

Appendices

43

Appendix A. Library Size Distribution Optimization

44

Appendix B. Limited Warranty

48

Preface

Regulatory Disclaimer

For Research Use Only.

Not for use in diagnostic procedures.

Contact Information

Technical Support

If you have questions, contact your local Roche Technical Support.



Go to sequencing.roche.com/support for contact information.

Manufacturing & Distribution

Manufacturer	Roche Diagnostics Cape Town, South Africa
Distribution	Roche Diagnostics GmbH Mannheim, Germany
Distribution in USA	Roche Diagnostics Corporation Indianapolis, IN USA

Conventions Used in This Manual

Symbols

Symbols	Description
	Important Note: Information critical to the success of the procedure or use of the product. Failure to follow these instructions could result in compromised data.
	Information Note: Designates a note that provides additional information concerning the current topic or procedure.

Text

Conventions	Description
Numbered listing	Indicates steps in a procedure that must be performed in the order listed.
<i>Italic type, blue</i>	Highlights a resource in a different area of this manual or on a web site.
<i>Italic type</i>	Identifies the external general resources or names.
Bold type	Identifies names of paragraphs, sections or emphasized words.



Chapter 1. Before You Begin


These Instructions for Use describe the process for depletion of ribosomal RNA (rRNA) and globin mRNA transcripts and streamlined RNA library construction of stranded RNA libraries for the Illumina sequencing ecosystem. Specifically, these Instructions for Use provide a protocol for the workflow outlined in [Figure 1](#) using the KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin.



The KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin Workflow provides:

- An easy to use, streamlined, and automation-friendly workflow with minimal resource requirements.
- Single vendor service when using the following accessory reagents:
 - KAPA Pure Beads
 - KAPA Unique Dual-Indexed (UDI) Adapters & KAPA Library Amplification Primer Mix or
 - KAPA Universal Adapter & KAPA UDI Primer Mixes 1 - 384

Figure 1 - KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin Quick Guide

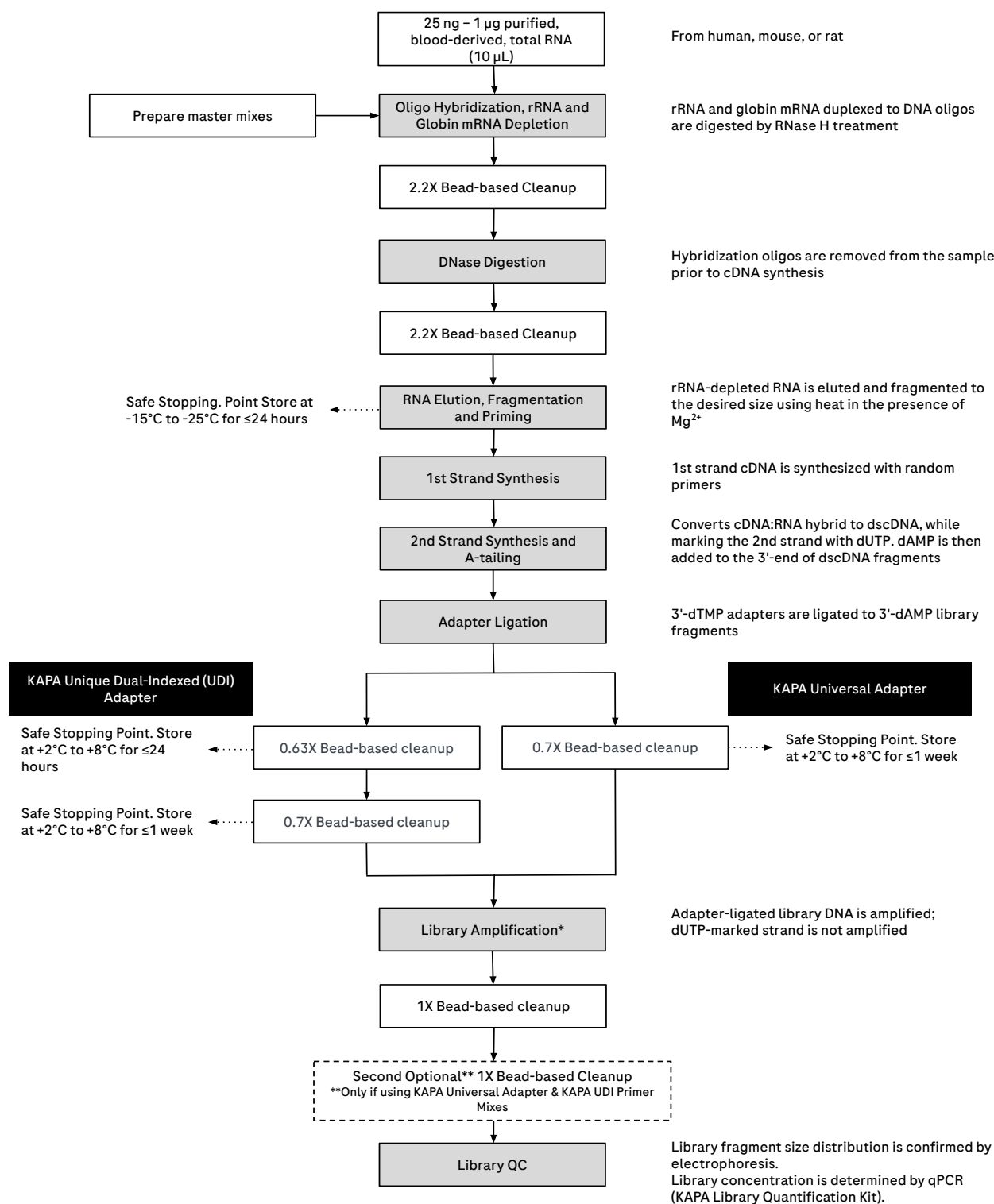

*Using KAPA HiFi HS ReadyMix and **either** KAPA Library Amplification Primer Mix if KAPA UDI Adapters were used **or** KAPA UDI Primers Mixes if KAPA Universal Adapter was used.

**A modified post-amplification cleanup is highly recommended when using KAPA Universal Adapter & KAPA UDI Primer Mixes with KAPA library preparation kits. Especially if libraries will be sequenced on an Illumina NovaSeq, NextSeq 1000/2000 or HiSeq X system (with patterned flow cells). This additional or modified cleanup will serve to remove any indexed primer carryover which may exacerbate index hopping.

To mix samples containing RNA, gently pipette the reaction mixture several times. Vortexing may fragment the RNA, resulting in lower quantity and a reduced library insert size.

Protocol Information & Safety

- Wear gloves and take precautions to avoid sample contamination.
- Perform all centrifugations at room temperature (+15°C to +25°C).
- Unless otherwise specified, all mixing steps are listed as 'mix thoroughly' and indicate that mixing should be performed by vortexing for at least 10 seconds.
- If liquid has collected in a tube's cap after mixing, gently tap or briefly centrifuge the sample to collect the liquid into the tube's bottom, ensuring that the mixture remains homogeneous before progressing to the next step.

References

- Thermocycler Manual
- Qubit Fluorometer Manual
- Qubit dsDNA HS Assay Kit Guide
- Qubit RNA HS Assay Kit Guide

Terminology

Sample Library: The initial shotgun library generated from RNA by fragmentation, 1st and 2nd strand synthesis, A-tailing and adapter ligation.

Amplified Sample Library: The sample library after amplification by PCR.

KAPA UDI Adapter: KAPA Unique Dual-Indexed Adapters are full-length, ready-to-use, QC-tested, dual-indexed Illumina adapters for ligation-based library construction.

KAPA UDI Primer Mixes: KAPA Unique Dual-Indexed Primer Mixes to be used in conjunction with KAPA Universal Adapters.

KAPA Universal Adapter: Truncated adapter containing a subset of sequencing motifs. Used in conjunction with KAPA UDI Primer Mixes.


Prepare the Following Equipment and Reagents

- Thermocyclers should be programmed with the following:
 - Oligo Hybridization and rRNA/globin mRNA Depletion program ([Chapter 3, Prepare the Sample Library, Step 1](#))
 - DNase Digestion program ([Chapter 3, Prepare the Sample Library, Step 3](#))
 - RNA Elution, Fragmentation and Priming program ([Chapter 3, Prepare the Sample Library, Step 5](#))
 - 1st Strand Synthesis program program ([Chapter 3, Prepare the Sample Library, Step 6](#))
 - 2nd Strand Synthesis and A-tailing program ([Chapter 3, Prepare the Sample Library, Step 7](#))
 - Adapter Ligation program ([Chapter 3, Prepare the Sample Library, Step 8](#))
 - Library Amplification program ([Chapter 4, Amplify the Sample Library, Step 2](#))

 It is recommended to use a thermocycler with a programmable heated lid. If not possible, please use the default settings.

- The following steps should be taken before beginning the workflow:

 Verify you are using the most up-to-date version of these Instructions for Use, go to sequencing.roche.com/support

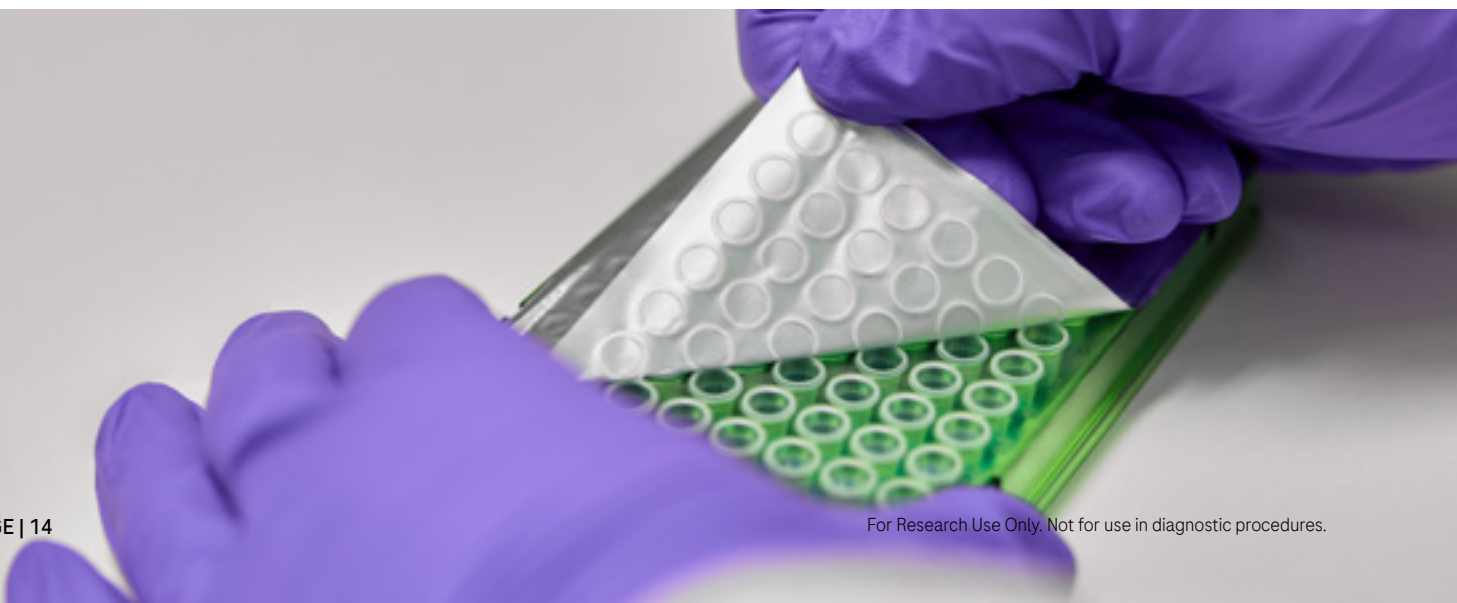

Required Equipment, Labware & Consumables

You assume full responsibility when using the equipment, labware, and consumables described below. This protocol is designed for use with the specified equipment, labware, and consumables.

Laboratory Equipment

Equipment	Supplier
DynaMag-96 Side Magnet	Thermo Fisher
Microcentrifuge (16,000 x g capability)	Multiple Vendors
Qubit Fluorometer	ThermoFisher
Electrophoretic device & associated assays and reagents	Multiple Vendors
Thermocycler	Multiple Vendors
Vortex mixer	Multiple Vendors
Plate Centrifuge (minimum 280 x g capability)	Multiple Vendors

Consumables Available from Roche

Description	Package Size	Material Number
KAPA Library Quantification Kit	Various options	Various material numbers
KAPA Pure Beads	5 mL	07983271001
	30 mL	07983280001
	60 mL	07983298001
KAPA PEG/NaCl Solution	20 mL	07961928001
KAPA Unique Dual-Indexed Adapter Kit	96 x 20 µL	08861919702
KAPA Library Amplification Primer Mix	250 reactions (1.25 mL)	07958994001
KAPA Library Amplification Primer Mix	384 reactions (1.92 mL)	09420410001
KAPA Library Amplification Primer Mix 96-well plate	96 x 5 µL	09420479001
KAPA Universal Adapter	96 reactions	09063781001
	384 reactions*	09063790001
KAPA UDI Primer Mixes 1 - 96	96 reactions	09134336001
KAPA UDI Primer Mixes 97 - 192	96 reactions	09329838001
KAPA UDI Primer Mixes 193 - 288	96 reactions	09329846001
KAPA UDI Primer Mixes 289 - 384	96 reactions	09329854001

* Virtual kits - consist of 4 x 96 reaction kits



Consumables Purchased from Other Vendors

Component	Supplier
10 mM Tris-HCl, pH 8.0 – 8.5	Multiple Vendors
Ethanol, 200 proof (absolute), for molecular biology	Multiple Vendors
Qubit dsDNA HS Assay Kit & Qubit RNA HS Assay Kit	ThermoFisher
Qubit Assay Tubes	ThermoFisher
Low binding Tubes: <ul style="list-style-type: none">▪ 0.2 mL PCR tubes▪ 1.5 mL microcentrifuge tubes (optional)	Multiple Vendors
Nuclease-free, PCR-grade water	Multiple Vendors
Agilent RNA 6000 Pico Kit	Agilent


Chapter 2.

Store and Prepare the Reagents




Step 1. Store the Reagent Kits

Reagent Kit	Storage Temperature
KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin	-15°C to -25°C
KAPA Pure Beads*	+2°C to +8°C
KAPA UDI Adapter Kit	-15°C to -25°C
KAPA Universal Adapter	-15°C to -25°C
KAPA UDI Primer Mixes or KAPA UDI Primer Mixes (resuspended)	+2°C to +8°C or -15°C to -25°C


 *The KAPA Pure Beads must not be frozen.

Step 2. Prepare KAPA Adapters and KAPA UDI Primer Mixes


 For multiplexing guidelines, please refer to either the *KAPA UDI Adapter* or *KAPA UDI Primer Mixes Instructions for Use* (available [online](#) on eLabDoc) as applicable.


Step 2a. Dilute the KAPA Unique Dual-Indexed Adapters (if required)


1. Retrieve the KAPA UDI Adapter plate from storage (-15°C to -25°C) and thaw at room temperature.
2. Centrifuge the KAPA UDI Adapter plate at room temperature (280 x g for at least 1 min) to ensure the liquid is collected at the bottom of the wells.

 Do not vortex the adapter plate as it could result in cross-contamination of the KAPA UDI Adapters. Pipette-mix individual adapters prior to use.

3. Before removing the foil cover, please ensure the plate is in the correct orientation. Please refer to the KAPA UDI Adapter Technical Data Sheet (Document number KR1736, available [online](#) on eLabDoc) for additional handling instructions.
4. Upon first use, carefully remove the foil cover of the plate to avoid cross contamination.

 Discard the original foil cover. Do not reuse.

 Maintain good lab practices when removing the foil seal. E.g., refrain from contact with the underside of the seal. If contact is made, change gloves before handling the adapter or sample plates.

 A new pipette tip must be used for each well to avoid cross contamination. If you are not using the entire contents of the KAPA UDI Adapter plate at this time, apply a new adhesive foil seal provided in the kit. Make sure that the foil is properly aligned and fully covers all 96 wells. Use a roller or other appropriate tool to ensure that the foil is evenly applied.

5. If applicable, dilute adapters to the required concentration, see Table 1 below, using the KAPA Adapter Dilution Buffer.

Step 2b. Dilute the KAPA Universal Adapters (if required)

1. Retrieve the KAPA Universal Adapters tube from storage (-15°C to -25°C) and thaw at room temperature or on ice.
2. Mix thoroughly before first use.
3. Centrifuge the KAPA Universal Adapter at room temperature (280 x g for at least 1 min) to ensure the liquid is collected at the bottom of the tube.
4. If applicable, dilute KAPA Universal Adapters to the required concentration, see Table 1 below, using 10 mM Tris-HCl pH 8.0 - 8.5.

Table 1. Recommended adapter concentrations for libraries constructed from different RNA sample quantities

Quantity of starting material	Adapter stock concentration *
25 – 499 ng	1.5 µM
500 – 1000 ng	7 µM

*Concentration required when using both KAPA UDI Adapter (full length) and KAPA Universal Adapter (truncated). Both KAPA UDI Adapter and KAPA Universal Adapter are provided at 15 µM.

Step 2c. Preparation of the KAPA UDI Primer Mixes

Before use of the KAPA UDI Primer Mixes, undertake the following steps to resuspend the primers:

1. Retrieve the KAPA UDI Primer Mixes plate from storage (+2°C to +8°C).
2. Centrifuge the KAPA UDI Primer Mixes plate at 280 x g for 1 minute to ensure the contents are at the bottom of the wells.
3. Before removing the foil cover, please ensure the plate is in the correct orientation before proceeding. In order to have well position A1 on the top left corner, the notched corner must be facing the user on the bottom left, as shown in [Figure 2](#).
4. Carefully remove the foil cover on the plate ensuring to avoid cross contamination. Discard the original foil cover.
5. Using a multichannel pipette, add 10 µL of Nuclease-free, PCR-grade water directly to the bottom of each well and discard tips after dispensing.



A new pipette tip must be used for each well to avoid cross contamination. Be sure to dispense water slowly to the bottom of each well to avoid liquid splash over to adjacent wells.

6. Ensure every well contains 10 µL of Nuclease-free, PCR-grade water and cover the plate with one of the adhesive foil seals provided in the kit.



Make sure the foil seal is properly aligned and fully covers all 96 wells. Failure to do so can lead to cross contamination of the KAPA UDI Primer Mixes.

7. Use a roller or appropriate tool to ensure the foil seal is evenly applied.
8. Centrifuge the plate at 280 x g for 30 seconds to ensure the dispensed 10 µL are at the bottom of the well.
9. Thoroughly vortex the plate ensuring all wells are mixed well.



Ensure wells at the corners of the plate are mixed well by vortexing the corners of the plate. Keep the plate upright.

10. Centrifuge the plate at 280 x g for 1 minute to ensure the contents are collected at the bottom of the wells.

11. The KAPA UDI Primer Mixes plate is now ready for use in the *Library Amplification* step
12. Store any unused but already resuspended KAPA UDI Primer Mixes at -15°C to -25°C. To avoid repeated freeze/thaw cycles you may transfer the resuspended primers to separate tubes or tube strips for storage

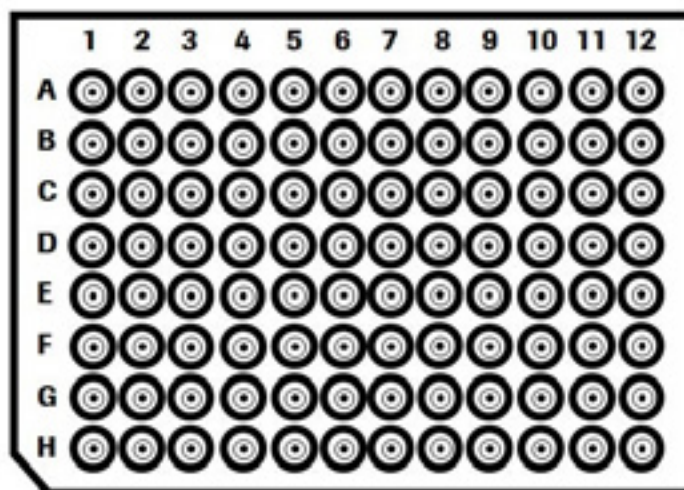


Figure 2: KAPA UDI Primer Mixes plate layout.

Considerations for Adapter Design and Concentration

- Adapter concentration affects ligation efficiency, as well as adapter and adapter-dimer carry-over during the post-ligation cleanup. The optimal adapter concentration for your workflow represents a compromise between the above factors and cost.
- When optimizing workflows for RNA inputs ≤ 25 ng, two or three adapter concentrations should be evaluated: try the recommended adapter concentration, as well as one or two additional concentrations in a range that is 2 – 10 times higher than the recommended concentration.
- Adapter quality has an impact on the effective concentration of adapter available for ligation. Always source the highest quality adapters from a reliable supplier, dilute and store adapters in a buffered solution with the requisite ionic strength, and avoid excessive freezing and thawing of adapter stock solutions.
- To accommodate different adapter concentrations within a batch of samples processed together, it is best to vary the concentrations of adapter stock solutions, and dispense a fixed volume (5 μ L) of each adapter. The alternative (using a single stock solution, and dispensing variable volumes of adapter into ligation reactions) is not recommended.

Chapter 3.

Prepare the Sample Library



Product Description

The KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin contains all of the buffers and enzymes required for depletion of ribosomal RNA (rRNA) and globin mRNA transcripts and the rapid construction of stranded RNA-Seq libraries from 25 ng – 1 µg of purified total RNA via the following steps:

1. depletion of rRNA and globin mRNA by hybridization of complementary DNA oligonucleotides, followed by treatment with RNase H and DNase to remove rRNA and globin mRNA duplexed to DNA and original DNA oligonucleotides, respectively;
2. fragmentation using heat and magnesium;
3. 1st strand cDNA synthesis using random priming;
4. combined 2nd strand synthesis and A-tailing, which converts the cDNA:RNA hybrid to double-stranded cDNA (dscDNA), incorporates dUTP into the second cDNA strand for stranded RNA sequencing, and adds dAMP to the 3' ends of the resulting dscDNA;
5. adapter ligation, where dsDNA adapters with 3' dTMP overhangs are ligated to library insert fragments; and
6. library amplification, to amplify library fragments carrying appropriate adapter sequences at both ends using high-fidelity, low-bias PCR. The strand marked with dUTP is not amplified, allowing strand-specific sequencing.

The kit provides KAPA Pure Beads for reaction cleanups, but does not include RNA or adapters. KAPA Adapters and KAPA UDI Primer Mixes are sold separately.

Reaction buffers are supplied in convenient formats comprising all of the required reaction components. This is designed to minimize the risk of RNase contamination, maximize homogeneous reaction composition and improve uniformity among replicate samples. Similarly, a single enzyme mixture is provided for each step of the library construction process, reducing the number of pipetting steps.

In order to maximize sequence coverage uniformity and to maintain relative transcript abundance, it is critical that library amplification bias be kept to a minimum. KAPA HiFi DNA Polymerase has been designed for low-bias, high-fidelity PCR, and is the polymerase of choice for NGS library amplification^{1,2,3,4}. The KAPA RNA HyperPrep Kit includes KAPA HiFi HotStart ReadyMix (2X) and KAPA Library Amplification Primer Mix (10X) for library amplification.

1. Oyola, S.O., et al., BMC Genomics 13, 1 (2012).
2. Quail, M.A., et al., Nature Methods 9, 10 – 11 (2012).
3. Quail, M.A., et al., BMC Genomics 13, 341 (2012).
4. Ross, M.G., et al., Genome Biology 14, R51 (2013).

The workflow requires the use of components from the following kits:

- KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin
- KAPA Unique Dual-Indexed Adapter Kit & KAPA Library Amplification Primers OR*
- KAPA Universal Adapter** & KAPA UDI Primer Mixes

Ensure that the following are available:

- Nuclease-free, PCR-grade water
- Freshly-prepared 80% ethanol
- 10 mM Tris-HCl, pH 8.0 - 8.5

*KAPA Unique Dual-Indexed (UDI) Adapters are recommended for use with the KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin. However, the kit is also compatible with KAPA Universal Adapter & KAPA UDI Primer Mixes as well as other full-length or truncated adapter designs from 3rd party suppliers.

**The KAPA Universal UMI Adapter (Catalog# 09329862001, and provided at 33 µM), may also be considered for this workflow. The KAPA Universal UMI Adapter with Unique Molecular Identifiers (UMI) allows for molecular barcoding within each sample, enabling proper molecule counting for low-frequency variant detection. However, the KAPA Universal UMI Adapter has only been validated within specific KAPA Target Enrichment workflows and may not confer the same computational benefits for non-target enrichment workflows (such as whole transcriptome sequencing). Contact your local Roche Technical Support at sequencing.roche.com/support for guidance and important considerations.

Sample Requirements

- The protocol has been validated for library construction from 25 ng – 1 µg of blood-derived, purified total RNA in ≤10 µL of RNase-free water. The quantity of rRNA and globin mRNA transcripts in a total RNA sample can vary significantly between samples. An input of 25 ng – 1 µg of total RNA is recommended to ensure that sufficient rRNA/globin mRNA-depleted total RNA is available for downstream library preparation.
- RNA in volumes >10 µL should be concentrated to 10 µL prior to use by either ethanol precipitation, bead purification (e.g., KAPA Pure Beads or RNAClean XP beads, Beckman Coulter), or column-based methods (e.g., RNeasy MinElute Cleanup Kit, QIAGEN). **Note that some loss of material is inevitable when using any of the above methods to concentrate RNA.**
- When concentrating RNA, elute in 12 µL of RNase-free water to ensure that 10 µL is available for use with this protocol.
- It is recommended to assess the quality and size distribution of the input RNA prior to rRNA/globin mRNA depletion by an electrophoretic method (e.g., Agilent Bioanalyzer RNA assay).
- In some cases, inhibitors may be present in a blood-derived total RNA sample that could limit the efficiency of the rRNA and globin mRNA depletion. If the presence of inhibitors is unknown or suspected, an optional, up-front 3X bead purification may be used (e.g., KAPA Pure Beads or RNAClean XP, Beckman Coulter) to remove any inhibitors. Please refer to the relevant Technical Data Sheet (Document number KR1245, available [online](#) on eLabDoc) for a detailed cleanup protocol.
 - Please be advised that the KAPA Pure Beads volume supplied in this kit may not be sufficient for the implementation of a 3X cleanup for all samples. If a 3X cleanup is deemed necessary, please purchase additional reagents.

RNA Handling

- RNases are ubiquitous and special care should be taken throughout the procedure to avoid RNase contamination.
- To avoid airborne RNase contamination, keep all reagents and RNA samples closed when not in use.
- Use a laminar flow hood if available, or prepare a sterile and RNase-free area. Clean the workspace, pipettes and other equipment with an RNase removal product (e.g., RNaseZAP, Ambion Inc.) according to manufacturer's recommendations.
- To avoid RNase contamination, always wear gloves when handling reagents, and use certified RNase-free plastic consumables. Change gloves after making contact with equipment or surfaces outside of the RNase-free working area.
- To mix samples containing RNA, gently pipette the reaction mixture several times. Vortexing may fragment the RNA, resulting in lower quantity and a reduced library insert size.
- To avoid degradation, minimize the number of freeze-thaw cycles and always store RNA in RNase-free water.

Important considerations for RNA Fragmentation

- RNA is fragmented by incubating at a high temperature in the presence of magnesium before carrying out 1st strand cDNA synthesis.
- Fragmentation conditions given in [Chapter 3](#) should be used as a guideline and may require adjustment based upon the quality and size distribution of the input RNA. It is recommended that a non-precious, representative sample of RNA be evaluated for the optimal fragmentation conditions.
- For fragmentation optimization beyond what is provided in [Chapter 3](#), please refer to [Appendix A: Library Size Distribution Optimization](#).

Step 1. Oligo Hybridization and rRNA/globin mRNA Depletion

This protocol requires 25 ng – 1 µg of purified, blood-derived, total RNA, in 10 µL of RNase-free water.



Ensure that the hybridization master mix and the depletion master mix are prepared and kept at room temperature before use.



If also processing non blood-derived samples, or if following the KAPA RiboErase (HMR) workflow, the volume of Globin Hybridization Oligos (HMR) must be replaced with equal volume RNase-free water. For more information, refer to the *KAPA RNA HyperPrep Kit with RiboErase (HMR) Instructions for Use* (available [online](#) on eLabdoc).

1. Pre-program a thermocycler as per the table below. Ensure that the thermocycler lid temperature is set to +105°C (default).

Step	Temperature	Duration
Hybridization	+95°C	2 min
Ramp down to 45°C at -0.1°C/s		
PAUSE	+45°C	∞
Depletion	+45°C	30 min
Hold	+4°C	∞

2. Prepare the required volume of hybridization master mix and depletion master mix as follows and keep at room temperature before use:

Oligo Hybridization Master Mix

Component	Volume Per Individual Sample
Globin Hybridization Oligos (HMR)*	1 µL
Hybridization Buffer	4 µL
Hybridization Oligos (HMR)	4 µL
RNase-free water	1 µL
Total hybridization master mix volume	10 µL

*If also processing non blood-derived samples, or if following the KAPA RiboErase (HMR) workflow, the volume of Globin Hybridization Oligos (HMR) must be replaced with an equal volume of RNase-free water (1 µL). For more information refer to the *KAPA RNA HyperPrep Kit with RiboErase (HMR) Instructions for Use* (available [online](#) on eLabdoc)


Depletion Master Mix

Component	Volume Per Individual Sample
Depletion Buffer	3 µL
RNase H	2 µL
Total depletion master mix volume	5 µL

3. Assemble the rRNA/globin mRNA hybridization reaction as per table below, and place samples in the pre-programmed thermocycler and execute the program:

Component	Volume Per Individual Sample
Total RNA in water	10 µL
Hybridization master mix (room temperature)	10 µL
Total volume	20 µL

- When the program reaches the pause step at +45°C, add the Depletion master mix to each 20 µL hybridization reaction and mix thoroughly by pipetting up and down multiple times.

 Ensure the depletion master mix containing RNase H is added while the samples are kept at +45°C in the thermocycler. Do NOT remove the samples from the thermocycler.

Component	Volume Per Individual Sample
Total RNA hybridized to oligos (at +45°C)	20 µL
Depletion master mix (room temperature)	5 µL
Total volume	25 µL


- Resume the cycling program to continue with the depletion step (+45°C for 30 minutes).
- Proceed immediately to [Step 2. rRNA/globin mRNA Depletion Cleanup](#)

Step 2. rRNA/globin mRNA Depletion Cleanup


- To each rRNA/globin mRNA depleted reaction, add 55 µL of room temperature KAPA Pure Beads that have been thoroughly resuspended as per the table below:

Component	Volume Per Individual Sample
rRNA/globin mRNA-depleted total RNA	25 µL
KAPA Pure Beads	55 µL
Total volume	80 µL

- Once added, pipette-mix thoroughly and centrifuge briefly to collect all droplets. Do NOT allow beads to pellet.

 It is important at this step to ensure that the solution is thoroughly mixed and appears homogeneous.

- Incubate the sample at room temperature for 5 min to bind RNA to the beads.
- Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
- Carefully remove and discard 75 µL of supernatant.
- Keeping the sample on the magnet, add 200 µL of 80% ethanol.
- Incubate the sample on the magnet at room temperature for ≥30 sec.
- Carefully remove and discard the ethanol.
- Keeping the sample on the magnet, add 200 µL of 80% ethanol.
- Incubate the sample on the magnet at room temperature for ≥30 sec.
- Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.
- Dry the beads at room temperature for 3 – 5 min, or until all of the ethanol has evaporated.

 Over-drying the beads may result in dramatic yield loss. Over-drying is indicated by a dry, cracked appearance of the bead pellet. The bead pellet should have a matte appearance when sufficiently dried.

Step 3. DNase Digestion

To remove the hybridization oligonucleotides from the rRNA/globin mRNA-depleted RNA, the sample is incubated with DNase.



Ensure that the DNA digestion master mix is prepared and kept at room temperature.



Do not vortex the DNase or DNA digestion master mix. Mix thoroughly by pipetting up and down multiple times.

1. Prepare the required volume of DNase digestion master mix as follows and keep at room temperature before use:

DNase Digestion Master Mix

Component	Volume Per Individual Sample
DNase Buffer	2.2 μ L
DNase	2 μ L
RNase-free water	17.8 μ L
Total DNase digestion master mix	22 μL

2. Assemble the DNase Digestion reactions as follows:

Component	Volume Per Individual Sample
Beads with rRNA/globin mRNA-depleted RNA	-
DNase digestion master mix	22 μ L
Total volume	22 μL

3. Once added, pipette-mix thoroughly and centrifuge briefly to collect all droplets. Do **NOT** allow beads to pellet.



It is important at this step to ensure that the solution is thoroughly mixed and appears homogeneous.

4. Incubate the sample at room temperature for 3 min to elute the RNA off the beads.
5. Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
6. Carefully transfer 20 μ L of supernatant into a new plate/tube(s). Discard the plate/tube(s) with beads.
7. Incubate the sample(s) with supernatant using the following protocol. Ensure that the thermocycler lid temperature is set to +50°C.

Step	Temperature	Duration
DNase digestion	+37°C	30 min
Hold	+4°C	∞

8. Proceed immediately to [Step 4. DNase Digestion Cleanup](#).

Step 4. DNase Digestion Cleanup

1. To each DNase depleted reaction, add 44 μL of room temperature KAPA Pure Beads that have been thoroughly resuspended as per the table below:

Component	Volume Per Individual Sample
DNase-treated RNA	20 μL
KAPA Pure Beads	44 μL
Total volume	64 μL

2. Once added, pipette-mix thoroughly and centrifuge briefly to collect all droplets. Do NOT allow beads to pellet.



It is important at this step to ensure that the solution is thoroughly mixed and appears homogeneous.

3. Incubate the sample at room temperature for 5 min to bind RNA to the beads.
4. Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
5. Carefully remove and discard 60 μL of supernatant.
6. Keeping the sample on the magnet, add 200 μL of 80% ethanol.
7. Incubate the sample on the magnet at room temperature for ≥ 30 sec.
8. Carefully remove and discard the ethanol.
9. Keeping the sample on the magnet, add 200 μL of 80% ethanol.
10. Incubate the sample on the magnet at room temperature for ≥ 30 sec.
11. Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.
12. Dry the beads at room temperature for 3 – 5 min, or until all of the ethanol has evaporated.



Over-drying the beads may result in dramatic yield loss. Over-drying is indicated by a dry, cracked appearance of the bead pellet. The bead pellet should have a matte appearance when sufficiently dried.

Step 5. RNA Elution, Fragmentation and Priming

RNA depleted of rRNA and globin mRNA is eluted from beads in Fragment, Prime and Elute Buffer (1X) and fragmented to the desired size by incubation at high temperature.

1. Prepare the required volume of Fragment, Prime and Elute Buffer (1X) at room temperature as per the table below:

Component	Volume Per Individual Sample
Fragment, Prime and Elute Buffer (2X)	11 μ L
RNase-free water	11 μ L
Total volume	22 μL

2. Thoroughly resuspend the beads with purified, DNase-treated RNA in 22 μ L of Fragment, Prime and Elute Buffer (1X) by gently pipetting up and down multiple times.
3. Incubate the sample at room temperature for 3 min to elute RNA off the beads.
4. Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
5. Carefully transfer 20 μ L of supernatant to a new plate/tube(s). Discard the plate/tube(s) with beads.



Safe stopping point - Samples can be stored at -15°C to -25°C for \leq 24 hrs. When ready, proceed to the next step.

6. Place the plate/tube(s) in a thermocycler and carry out the fragmentation and priming program as per the table below:



Ensure that the thermocycler lid temperature is set to +105°C (default).

Desired mean library insert size (bp)	Fragmentation
100 - 200	8 min at +94°C
200 - 300	6 min at +94°C
300 - 400	6 min at +85°C

7. Place the plate/tube(s) on ice and proceed immediately to [Step 6. 1st Strand Synthesis](#).

Step 6. 1st Strand Synthesis

1. On ice, assemble the 1st strand synthesis master mix as per the table below:

1st Strand Synthesis Master Mix

Component	Volume Per Individual Sample <i>inc. 20% excess</i>
1st Strand Synthesis Buffer	11 μ L
KAPA Script	1 μ L
Total 1st strand synthesis master mix volume	12 μL

2. Keeping the plate/tube(s) on ice, combine the input RNA (fragmented & primed) and 1st strand synthesis master mix as per the table below:

Component	Volume Per Individual Sample
Fragmented, primed RNA	20 μ L
1st strand synthesis master mix	10 μ L
Total volume	30 μL

3. Keeping the plate/tube(s) on ice, mix thoroughly by gently pipetting the reaction up and down several times.



Ensure that the thermocycler lid temperature is set to +85°C.

4. Incubate the plate/tube(s) using the following protocol:

Step	Temperature	Duration
Primer extension	+25°C	10 min
1st strand synthesis	+42°C	15 min
Enzyme inactivation	+70°C	15 min
Hold	+4°C	∞

5. Place the plate/tube(s) on ice and proceed immediately to [Step 7. 2nd Strand Synthesis & A-tailing](#).

Step 7. 2nd Strand Synthesis & A-tailing

1. On ice, assemble the 2nd strand synthesis and A-tailing master mix as per table below:

2nd Strand Synthesis & A-tailing Master Mix

Component	Volume Per Individual Sample <i>inc. 10% excess</i>
2nd Strand Marking Buffer	31 μ L
2nd Strand Synthesis & A-Tailing Enzyme Mix	2 μ L
Total 2nd strand synthesis and A-tailing master mix volume	33 μL

2. Keeping the plate/tube(s) on ice, combine the 1st strand synthesis product and 2nd strand synthesis and A-tailing master mix as per the table below:

Component	Volume Per Individual Sample
1st strand synthesis product	30 μ L
2nd strand synthesis & A-tailing master mix	30 μ L
Total volume	60 μL

3. Keeping the plate/tube(s) on ice, mix thoroughly by gently pipetting the reaction up and down several times.




Ensure that the thermocycler lid temperature is set to +85°C.


4. Incubate the plate/tube(s) using the following protocol:


Step	Temperature	Duration
2nd strand synthesis	+16°C	30 min
A-tailing	+62°C	10 min
Hold	+4°C	∞

5. Place the plate/tube(s) on ice and proceed immediately to [Step 8. Adapter Ligation](#).

Step 8. Adapter Ligation

 KAPA UDI Adapter and KAPA Universal Adapter, unless specified, will hereafter be referred to as KAPA Adapters.

 KAPA Adapters must be added to each tube/well individually prior to addition of the Ligation reagents.

 If using KAPA Universal Adapter, sample indexes are incorporated during the Library Amplification step ([Chapter 4, Amplify the Sample Library](#)). Precautions should be taken to avoid sample cross contamination.

1. On ice, assemble each Adapter Ligation reaction as per table below:

Adapter Ligation Master Mix

Component	Volume Per Individual Sample <i>inc. 10% excess</i>
Ligation Buffer	40 μ L
DNA Ligase	10 μ L
Total Adapter Ligation master mix volume	50 μL


2. Keeping the plate/tube(s) on ice, combine the 2nd strand synthesis product and diluted KAPA Adapter stock, then add the Adapter Ligation master mix as per the table below:

Component	Volume Per Individual Sample
2nd strand synthesis product	60 μ L
Adapter ligation master mix	45 μ L
Diluted adapter stock (Chapter 2)	5 μ L
Total volume	110 μL

3. Keeping the plate/tube(s) on ice, mix thoroughly by pipetting the reaction up and down several times.

 Ensure that the thermocycler lid temperature is set to +50°C.

4. Incubate the plate/tube(s) at +20°C for 15 min.
5. Proceed immediately to the next step.

 To achieve higher ligation efficiency, particularly for difficult samples, consider increasing the ligation time to a maximum of 4 hrs at +20°C or overnight at +2°C to +8°C. Please note that longer ligation times may lead to increased levels of adapter-dimer. Adapter concentrations may have to be optimized if ligation times are extended significantly.


6. Proceed immediately to [Step 9. Purify the Sample Library using KAPA Pure Beads](#).

Step 9. Purify the Sample Library using KAPA Pure Beads

 **If KAPA UDI Adapters were used for library construction – proceed to Step 9a or,**

 **If KAPA Universal Adapter were used for library construction – proceed to Step 9b**


Step 9a. Purify the Sample Library constructed using KAPA UDI Adapter (first cleanup)

 Only perform this 0.63X cleanup if full length, KAPA UDI Adapters were used. Proceed to [Step 9b](#) if the KAPA Universal Adapter was used in order to perform a single 0.7X post-ligation cleanup.


- To each ligation reaction, add 70 μL of room temperature KAPA Pure Beads that have been thoroughly resuspended as per the table below:

Component	Volume Per Individual Sample
Adapter-ligated DNA	110 μL
KAPA Pure Beads	70 μL
Total volume	180 μL


- Once added, mix thoroughly and centrifuge briefly to collect all droplets. Do **NOT** allow beads to pellet.

 It is important at this step to ensure that the solution is thoroughly mixed and appears homogeneous. Insufficient mixing may compromise recovery and result in size selection.


- Incubate the sample at room temperature for 5 min to bind DNA to the beads.
- Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
- Carefully remove and discard the supernatant.
- Keeping the sample on the magnet, add 200 μL of 80% ethanol.
- Incubate the sample on the magnet at room temperature for ≥ 30 sec.
- Carefully remove and discard the ethanol.
- Keeping the sample on the magnet, add 200 μL of 80% ethanol.
- Incubate the sample on the magnet at room temperature for ≥ 30 sec.
- Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.
- Dry the beads at room temperature for 3 – 5 min, or until all of the ethanol has evaporated.

 Over-drying the beads may result in dramatic yield loss. Over-drying is indicated by a dry, cracked appearance of the bead pellet. The bead pellet should have a matte appearance when sufficiently dried.

- Remove the sample from the magnet.
- Thoroughly resuspend the beads in 50 μL of elution buffer (10 mM Tris-HCl, pH 8.0 – 8.5) to proceed with a 2nd post-ligation cleanup using PEG/NaCl ([Step 10. Purify the Sample Library using PEG/NaCl](#)).
- Incubate the sample at room temperature for 2 minutes to allow the sample library to elute off the beads.

 **Safe Stopping Point** - If necessary this is a safe stopping point. The solution with resuspended beads can be stored at +2°C to +8°C for ≤ 24 hrs. **Do not freeze the beads, as this can result in dramatic loss of DNA.** When ready, proceed to [Step 10. Purify the Sample Library using PEG/NaCl](#).


Step 9b. Purify the Sample Library constructed using KAPA Universal Adapter

 Only perform this 0.7X cleanup if KAPA Universal Adapters were used. Refer to [Step 9a](#) if the KAPA UDI Adapters were used in order to perform a 0.63X post-ligation cleanup.


- To each ligation reaction, add 77 μL of room temperature KAPA Pure Beads that have been thoroughly resuspended as per the table below:

Component	Volume Per Individual Sample
Adapter-ligated DNA	110 μL
KAPA Pure Beads	77 μL
Total volume	187 μL

- Once added, mix thoroughly and centrifuge briefly to collect all droplets. Do **NOT** allow beads to pellet.

 It is important at this step to ensure that the solution is thoroughly mixed and appears homogeneous. Insufficient mixing may compromise recovery and result in size selection.

- Incubate the sample at room temperature for 5 min to bind DNA to the beads.
- Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
- Carefully remove and discard the supernatant.
- Keeping the sample on the magnet, add 200 μL of 80% ethanol.
- Incubate the sample on the magnet at room temperature for ≥ 30 sec.
- Carefully remove and discard the ethanol.
- Keeping the sample on the magnet, add 200 μL of 80% ethanol.
- Incubate the sample on the magnet at room temperature for ≥ 30 sec.
- Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.
- Dry the beads at room temperature for 3 – 5 min, or until all of the ethanol has evaporated.


 Over-drying the beads may result in dramatic yield loss. Over-drying is indicated by a dry, cracked appearance of the bead pellet. The bead pellet should have a matte appearance when sufficiently dried.

- Remove the sample from the magnet.
- Thoroughly resuspend the beads in 22 μL of 10 mM Tris-HCl, pH 8.0 – 8.5.
- Incubate the sample at room temperature for 2 minutes to allow the sample library to elute off the beads.
- Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
- Transfer 20 μL of the clear supernatant to a new plate/tube(s) and proceed to [Library Amplification \(Chapter 4, Step 1\)](#).



Safe Stopping Point - If necessary, this is a safe stopping point. Purified, adapter-ligated library may be stored at +2°C to +8°C for 1 – 2 weeks or at -15°C to -25°C for ≤ 1 month before amplification. To avoid degradation, always store DNA in a buffered solution (10 mM Tris-HCl, pH 8.0 – 8.5) when possible, and minimize the number of freeze-thaw cycles.


Step 10. Purify the Sample Library using PEG/NaCl (second cleanup)

 This is a **mandatory second cleanup** if KAPA UDI Adapters were used during library construction. Only perform this 0.7X cleanup if KAPA UDI Adapters were used. **SKIP** this step if KAPA Universal Adapters were used.

1. To each sample, add 35 μL of room temperature PEG/NaCl solution that have been thoroughly resuspended as per the table below:

Component	Volume Per Individual Sample
Beads with purified, adapter-ligated DNA (from Step 9a)	50 μL
PEG/NaCl Solution	35 μL
Total volume	85 μL

2. Once added, mix thoroughly and centrifuge briefly to collect all droplets. Do NOT allow beads to pellet.
3. Incubate the sample at room temperature for 5 min to bind DNA to the beads.
4. Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
5. Carefully remove and discard 80 μL of supernatant.
6. Keeping the sample on the magnet, add 200 μL of 80% ethanol.
7. Incubate the sample on the magnet at room temperature for ≥ 30 sec.
8. Carefully remove and discard the ethanol.
9. Keeping the sample on the magnet, add 200 μL of 80% ethanol.
10. Incubate the sample on the magnet at room temperature for ≥ 30 sec.
11. Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.
12. Dry the beads at room temperature for 3 – 5 min, or until all of the ethanol has evaporated.

 Over-drying the beads may result in dramatic yield loss. Over-drying is indicated by a dry, cracked appearance of the bead pellet. The bead pellet should have a matte appearance when sufficiently dried.

13. Remove the sample from the magnet.
14. Thoroughly resuspend the beads in 22 μL of 10 mM Tris-HCl (pH 8.0 – 8.5).
15. Incubate the sample at room temperature for 2 min to allow the sample library to elute off the beads.
16. Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
17. Transfer 20 μL of the clear supernatant to a new plate/tube(s) and proceed to [Library Amplification \(Chapter 4, Step 1\)](#).




Safe Stopping Point - If necessary, this is a safe stopping point. Purified, adapter-ligated library may be stored at +2°C to +8°C for 1 – 2 weeks or at -15°C to -25°C for ≤ 1 month before amplification. To avoid degradation, always store DNA in a buffered solution (10 mM Tris-HCl, pH 8.0 – 8.5) when possible, and minimize the number of freeze-thaw cycles.

Chapter 4.

Amplify the Sample Library



This chapter describes how to amplify the adapter-ligated library, using KAPA HiFi HotStart ReadyMix and appropriate compatible indexes or primer sets.

 If KAPA UDI Adapters were used for adapter ligation, ensure KAPA Library Amplification Primer Mix (provided with the kit) is used for library amplification.

 If KAPA Universal Adapters were used for adapter ligation, ensure a unique KAPA UDI Primer Mix is added to each sample library.







Ensure that the following is available:

- Freshly-prepared 80% ethanol
- 10 mM Tris-HCl, pH 8.0 – 8.5
- Nuclease-free, PCR-grade water

Library Amplification Considerations

- KAPA HiFi HotStart, the enzyme provided in the KAPA HiFi HotStart ReadyMix, is an antibody-based hotstart formulation of KAPA HiFi DNA Polymerase—a B-family DNA polymerase engineered for increased processivity and high fidelity. KAPA HiFi HotStart has 5'→3' polymerase and 3'→5' exonuclease (proofreading) activity, but no 5'→3' exonuclease activity. The error rate of KAPA HiFi HotStart is 2.8×10^{-7} errors/base, equivalent to 1 error in 3.5×10^6 nucleotides incorporated.
- The KAPA Library Amplification Primer Mix (10X) is designed to eliminate or delay primer depletion during library amplification reactions performed with KAPA HiFi HotStart ReadyMix. The primer mix is suitable for the amplification of all Illumina libraries flanked by the P5 and P7 flow cell sequences. Primers are supplied at a 10X concentration of 20 μ M each, and have been formulated as described below. User-supplied primer mixes may be used in combination with incomplete or custom adapters. Please contact Technical Support at sequencing.roche.com/support for guidelines on the formulation of user-supplied library amplification primers.
- To achieve the highest amplification efficiency and avoid primer depletion, it is critical to use an optimal concentration of high quality primers. Primers should be used at a final concentration of 0.5 – 4 μ M each.
- Library amplification primers should be HPLC-purified and modified to include a phosphorothioate bond at the 3'-terminal of each primer (to prevent degradation by the strong proofreading activity of KAPA HiFi HotStart DNA Polymerase). Always store and dilute primers in a buffered solution (e.g., 10 mM Tris-HCl, pH 8.0 – 8.5), and limit the number of freeze-thaw cycles. To achieve the latter, store primers at +2°C to +8°C for short-term use, or as single-use aliquots at -15°C to -25°C.
- In library amplification reactions (set up according to the recommended protocol), primers are typically depleted before dNTPs. When DNA synthesis can no longer take place due to substrate depletion, subsequent rounds of DNA denaturation and annealing result in the separation of complementary DNA strands, followed by imperfect annealing to non-complementary partners. This presumably results in the formation of so-called “daisy chains” or “tangled knots”, comprising large assemblies of improperly annealed, partially double-stranded, heteroduplex DNA. These species migrate slower and are observed as secondary, higher molecular weight peaks during the electrophoretic analysis of amplified libraries. However, they typically comprise library molecules of the desired length, which are individualized during denaturation prior to cluster amplification or probe hybridization. Since these heteroduplexes contain significant portions of single-stranded DNA, over-amplification leads to the under-quantification of library molecules with assays employing dsDNA-binding dyes. qPCR-based library quantifications methods, such as the KAPA Library Quantification assay, quantify DNA by denaturation and amplification, thereby providing an accurate measure of the amount of adapter-ligated molecules in a library—even if the library was over-amplified. Excessive library amplification can result in other unwanted artifacts such as amplification bias, PCR duplicates, chimeric library inserts and nucleotide substitutions. The extent of library amplification should therefore be limited as much as possible, while ensuring that sufficient material is generated for QC and downstream processing.
- If cycled to completion (not recommended), one 50 μ L library amplification can produce 8 – 10 μ g of amplified library. To minimize over-amplification and its associated, undesired artifacts, the number of amplification cycles should be tailored to produce the optimal amount of final library required for downstream processes. This is typically in the range of 250 ng – 1.5 μ g.

Step 1. Prepare the Library Amplification Reaction

-  Make sure the KAPA Pure Beads are removed from storage at least 30 minutes prior to starting the workflow to ensure they are equilibrated to room temperature. For best performance, store the beads protected from light when not in use.
-  KAPA HiFi HotStart ReadyMix (2X) contains isostabilizers and may not freeze completely, even when stored at -15°C to -25°C. Nevertheless, always ensure that the ReadyMix is **fully thawed** and **thoroughly mixed** before use.
-  If applicable, retrieve and thaw the KAPA UDI Primer Mixes plate prepared in [Chapter 2, Step 2c](#).
Centrifuge the KAPA UDI Primer Mixes plate at 280 x g for 30 seconds to collect the contents to the bottom of the wells and peel off or pierce the foil seal for the appropriate number of wells needed. If only using a subset of the KAPA UDI Primer Mixes from the original plate, remove and discard residual primers from the well and apply a new adhesive foil seal provided in the kit.
-  If piercing the foil seal, avoid cross contamination by using a new pipette tip for every well.
-  Proper re-sealing and storage of the KAPA UDI Primer Mixes plate is necessary for unused primer mixes utilization at a later date.
-  Keep all components on ice as long as possible during handling.

1. Assemble each library amplification reaction as per table below:

Component	Volume per Individual Library
KAPA HiFi HotStart ReadyMix (2X)	25 µL
KAPA Library Amplification Primer Mix* OR KAPA UDI Primer Mix**	5 µL
Adapter-ligated library	20 µL
Total volume	50 µL

* If KAPA UDI Adapters were used for adapter ligation, ensure KAPA Library Amplification Primer Mix is used for library amplification. The KAPA HiFi HotStart ReadyMix and KAPA Library Amplification Primer Mix should preferably be premixed and added in a single pipetting step.

** If KAPA Universal Adapters were used for adapter ligation ensure a unique KAPA UDI Primer Mix is added to each sample library.

2. Mix thoroughly and centrifuge briefly. Immediately proceed to the next step.

Step 2. Perform the Library Amplification

1. Place the sample in the thermocycler and amplify the Sample library using the following Library Amplification program with the lid temperature set to +105°C:

Step	Temperature	Duration	Cycles
Initial denaturation	+98°C	45 sec	1
Denaturation	+98°C	15 sec	Variable, see Table 2 below
Annealing	+60°C	30 sec	
Extension	+72°C	30 sec	
Final extension	+72°C	1 min	1
Hold	+4°C	∞	1

Table 2: Recommended cycle numbers

Quantity of starting material	Number of cycles*
25 - 100 ng	11 - 15
101 - 250 ng	9 - 12
251 - 500 ng	7 - 10
501 - 1000 ng	6 - 8

Note: The recommended number of cycles should be used as a guide for library amplification. Cycle numbers may require adjustment depending on RNA input quality, library amplification efficiency, presence of adapter-dimer, and the desired yield post amplification. Quantification of material after the second post-ligation cleanup using a qPCR assay, such as the KAPA Library Quantification Kit, can help to determine the number of amplification cycles required for a specific sample type or application.

*The recommended number of cycles is irrespective of which adapter kit was used during library construction (KAPA UDI Adapter or KAPA Universal Adapter).


2. Proceed immediately to the next step.

Step 3. Purify the Amplified Sample Library using KAPA Pure Beads

 **If KAPA UDI Adapters and KAPA Library Amplification Primer Mix were used for library construction – proceed to Step 3a or,**

 **If KAPA Universal Adapter and KAPA UDI Primer Mixes were used for library construction – proceed to Step 3b.**

Step 3a. Purify the Amplified Sample Library constructed using KAPA UDI Adapter & KAPA Library Amplification Primer Mix

1. Add 50 μL of room temperature, thoroughly resuspended, KAPA Pure Beads to each amplified sample library.
2. Mix the amplified sample library and KAPA Pure Beads thoroughly and centrifuge briefly to collect all droplets. Do NOT allow beads to pellet.
3. Incubate the sample at room temperature for 5 minutes to allow the amplified sample library to bind to the beads.
4. Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
5. Carefully remove and discard the supernatant.
6. Keeping the sample on the magnet, add 200 μL of freshly-prepared 80% ethanol.
7. Incubate the sample at room temperature for ≥ 30 seconds.
8. Carefully remove and discard the ethanol.
9. Keeping the sample on the magnet, add 200 μL of freshly-prepared 80% ethanol.
10. Incubate the sample at room temperature for ≥ 30 seconds.
11. Carefully remove and discard the ethanol. Remove residual ethanol without disturbing the beads.
12. Allow the beads to dry at room temperature, sufficiently for all of the ethanol to evaporate.
 Over-drying the beads may result in dramatic yield loss. Over-drying is indicated by a dry, cracked appearance of the bead pellet. The bead pellet should have a matte appearance when sufficiently dried.
13. Remove the sample from the magnet.
14. Thoroughly resuspend the beads in 22 μL (or appropriate volume) of 10 mM Tris-HCl, pH 8.0 – 8.5.
15. Incubate the sample at room temperature for 2 minutes to allow the amplified sample library to elute off the beads.
16. Place the sample on the magnet to capture the beads. Incubate until the liquid is clear.
17. Transfer the clear supernatant to a new tube(s)/well and proceed with library QC ([Chapter 5. Quality Control](#)), target capture or sequencing, as appropriate.
18. Purified, amplified sample libraries can be stored at +2°C to +8°C for 1 - 2 weeks or at -15°C to -25°C for up to 3 months.

Step 3b. Purify the Amplified Sample Library constructed using KAPA Universal Adapter & KAPA UDI Primer Mixes



A modified post-amplification cleanup is highly recommended when using KAPA Universal Adapter & KAPA UDI Primer Mixes with KAPA library preparation kits. Especially if libraries will be sequencing on an Illumina NovaSeq or HiSeq X system (with patterned flow cells). This additional or modified cleanup will serve to remove any indexed primer carryover which may exacerbate index hopping.



The KAPA Pure Beads provided with the 24 rxn and 96 rxn pack size, may not be sufficient if processing 24 or 96 samples with the KAPA Universal Adapter & KAPA UDI Primer Mixes. Additional KAPA Pure Beads is sold separately. Alternatively, the unused PEG/NaCl provided with the kit may be repurposed for a “with-bead” protocol. Please contact [Technical Support](#) for guidance.

1. Add 50 µL of room temperature, thoroughly resuspended, KAPA Pure Beads to each amplified sample library.
2. Mix the amplified sample library and KAPA Pure Beads thoroughly and centrifuge briefly to collect all droplets.
3. Incubate the sample at room temperature for 5 minutes to allow the amplified sample library to bind to the beads.
4. Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
5. Carefully remove and discard the supernatant.
6. Remove the tubes from the magnet, and resuspend the beads in 50 µL of nuclease-free PCR-grade water or 10 mM Tris-HCl, pH 8.0 - 8.5. Incubate for 2 minutes to allow the sample to elute off from the beads.
7. Add 50 µL of KAPA Pure Beads to each sample.
8. Mix thoroughly by pipetting or vortexing, and centrifuge briefly.
9. Incubate the sample at room temperature for 5 minutes to allow the amplified sample library to bind to the beads.
10. Place the sample on a magnet to capture the beads. Incubate until the liquid is clear.
11. Carefully remove and discard the supernatant.
12. Keeping the sample on the magnet, add 200 µL of freshly-prepared 80% ethanol.
13. Incubate the sample at room temperature for ≥30 seconds.
14. Carefully remove and discard the ethanol.
15. Keeping the sample on the magnet, add 200 µL of freshly-prepared 80% ethanol.
16. Incubate the sample at room temperature for ≥30 seconds.
17. Carefully remove and discard the ethanol. Remove residual ethanol without disturbing the beads.
18. Allow the beads to dry at room temperature, sufficiently for all of the ethanol to evaporate.



Over-drying the beads may result in dramatic yield loss. Over-drying is indicated by a dry, cracked appearance of the bead pellet. The bead pellet should have a matte appearance when sufficiently dried.

19. Remove the sample from the magnet.
20. Thoroughly resuspend the beads in 22 µL (or appropriate volume) of 10 mM Tris-HCl, pH 8.0 - 8.5.
21. Incubate the sample at room temperature for 2 minutes to allow the sample library to elute off the beads.
22. Place the sample on the magnet to capture the beads. Incubate until the liquid is clear.
23. Transfer the clear supernatant to a new tube(s)/well and proceed with library QC ([Chapter 5. Quality Control](#)), target capture or sequencing, as appropriate.
24. Purified, amplified sample libraries can be stored at +2°C to +8°C for 1 - 2 weeks or at -15°C to -25°C for up to 3 months.

Chapter 5.

Quality Control





This chapter describes how to determine the concentration and size distribution of the sample library. Your specific library construction workflow should be tailored and optimized to yield a sufficient amount of adapter-ligated molecules of the desired size distribution for the next step in the process (e.g., target capture or sequencing), as well as for library QC and archiving purposes.

Quantification

KAPA Library Quantification Kits are recommended for qPCR-based quantification of libraries generated using the KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin workflow. Libraries may also be quantified using a Qubit Fluorometer and the Qubit dsDNA HS Assay Kit.

Sizing

The size distribution of the final libraries should be confirmed with an electrophoretic method. A LabChip GX, GXII, or GX Touch (PerkinElmer), Bioanalyzer, TapeStation (Agilent Technologies), Fragment Analyzer (Advanced Analytical) system or similar instruments are recommended over conventional gels. Libraries generated using the KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin may require dilution prior to electrophoretic assessment. Refer to the respective instrument and assay user manuals for guidance on how to dilute your sample library so as not to exceed the detection limit of the specific assay.

Appendices

Appendix A. Library Size Distribution Optimization

The KAPA RNA HyperPrep Kit with RiboErase (HMR) Globin offers a tunable RNA fragmentation module in which RNA is fragmented at a high temperature in the presence of magnesium. Final library size distributions can be optimized for specific sample types and applications by varying both incubation time and temperature. Generally:

- higher temperatures and/or longer incubation times result in shorter, narrower distributions; and
- lower temperatures and/or shorter incubation times result in longer, broader distributions.

Fragmentation times may require adjustment based upon the quality of the input RNA. For intact RNA, such as that extracted from fresh/frozen tissue, longer fragmentation is required at higher temperatures. For degraded or fragmented RNA (e.g., from older samples), a lower temperature and/or shorter time should be used.

Intact Total RNA Inputs

When using high-quality, intact total RNA, the following recommendations can be used as a starting point for the optimization of final library distributions beyond what is provided within [Chapter 3](#). It is recommended that a non-precious, representative RNA sample be used for this optimization. For the following figures and tables, final libraries were generated using 100 ng of high-quality Universal Human Reference (UHR) RNA.

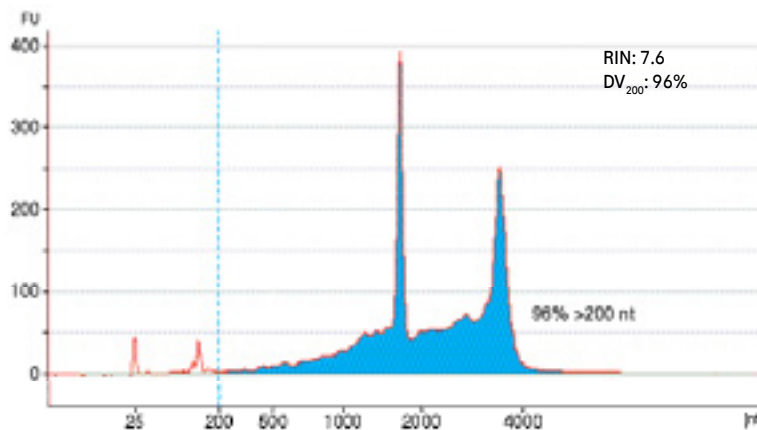


Figure 3. High-quality UHR total RNA electropherogram prior to library construction
RIN score was 7.6, with 96% of the RNA measuring >200 nucleotides, (i.e. the DV₂₀₀ score) as assessed via an Agilent RNA 6000 Pico Kit.

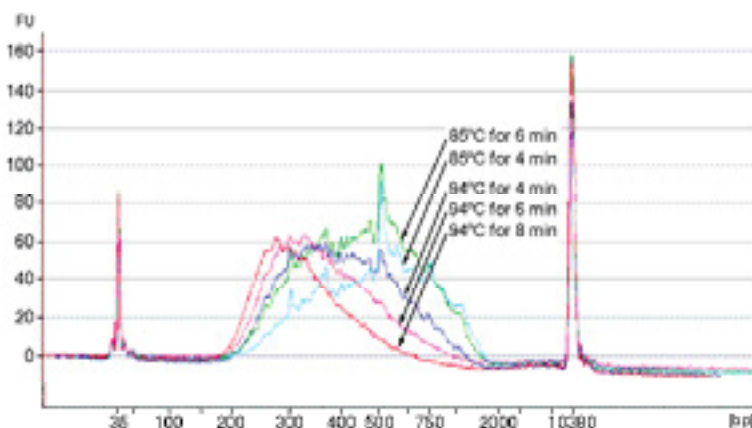
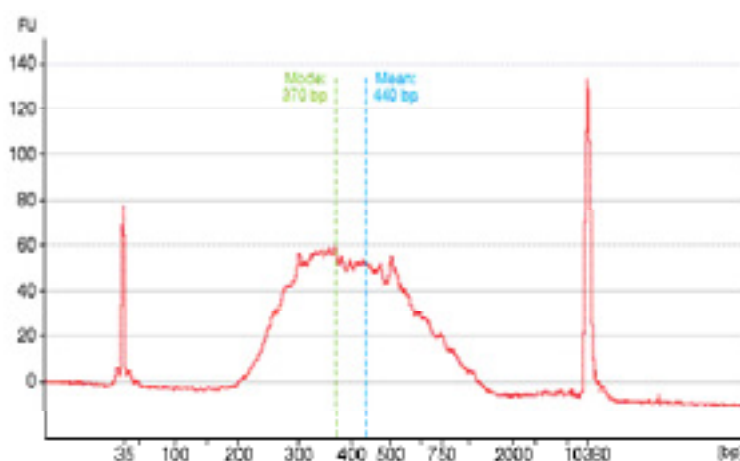


Figure 4. Final library distributions demonstrating fragmentation tunability.

Higher temperatures and longer incubation times resulted in shorter, narrower final library distributions. Libraries were constructed using 100 ng of high-quality UHR RNA and various fragmentation conditions. Note that results may differ with other sample sources. Electropherograms were generated with an Agilent High Sensitivity DNA Kit.

Table 3. Approximate mean and mode final library sizes (bp) for each fragmentation condition assessed


Fragmentation	Final library size (bp)	
	Mean	Mode
+94°C for 8 min	~350	~280
+94°C for 6 min	~380	~320
+94°C for 4 min	~440	~370
+94°C for 2 min	~550	~500
+85°C for 6 min	~510	~510
+85°C for 4 min	~570	~510

**Figure 5. Visual depiction of the mean and mode distribution metrics.**

For a final library generated using 100 ng UHR fragmented at +94°C for 4 minutes, the mode is the highest peak in the library (~370 bp), while the mean is the numerical average across all molecular lengths in the library (~440 bp). In this example, the mean of the library is calculated across the range of 190 to 1600 bp. The higher molecular weight shoulder of the distribution results in the mean being larger than the mode.

Partially Degraded Total RNA Inputs

When working with partially degraded (PD) inputs, the following recommendations may serve as a starting point for the optimization of final library distributions beyond what is provided within [Chapter 3](#). It is recommended that a non-precious, representative RNA sample be used for optimization. For the following figures and tables, chemically degraded samples of UHR were used as substitutes for real-world partially degraded samples of varying qualities. This RNA was ribosomally-depleted using the KAPA RiboErase (HMR) RNA enrichment module, processed with the KAPA RNA HyperPrep workflow, and subjected to a fragmentation condition titration.

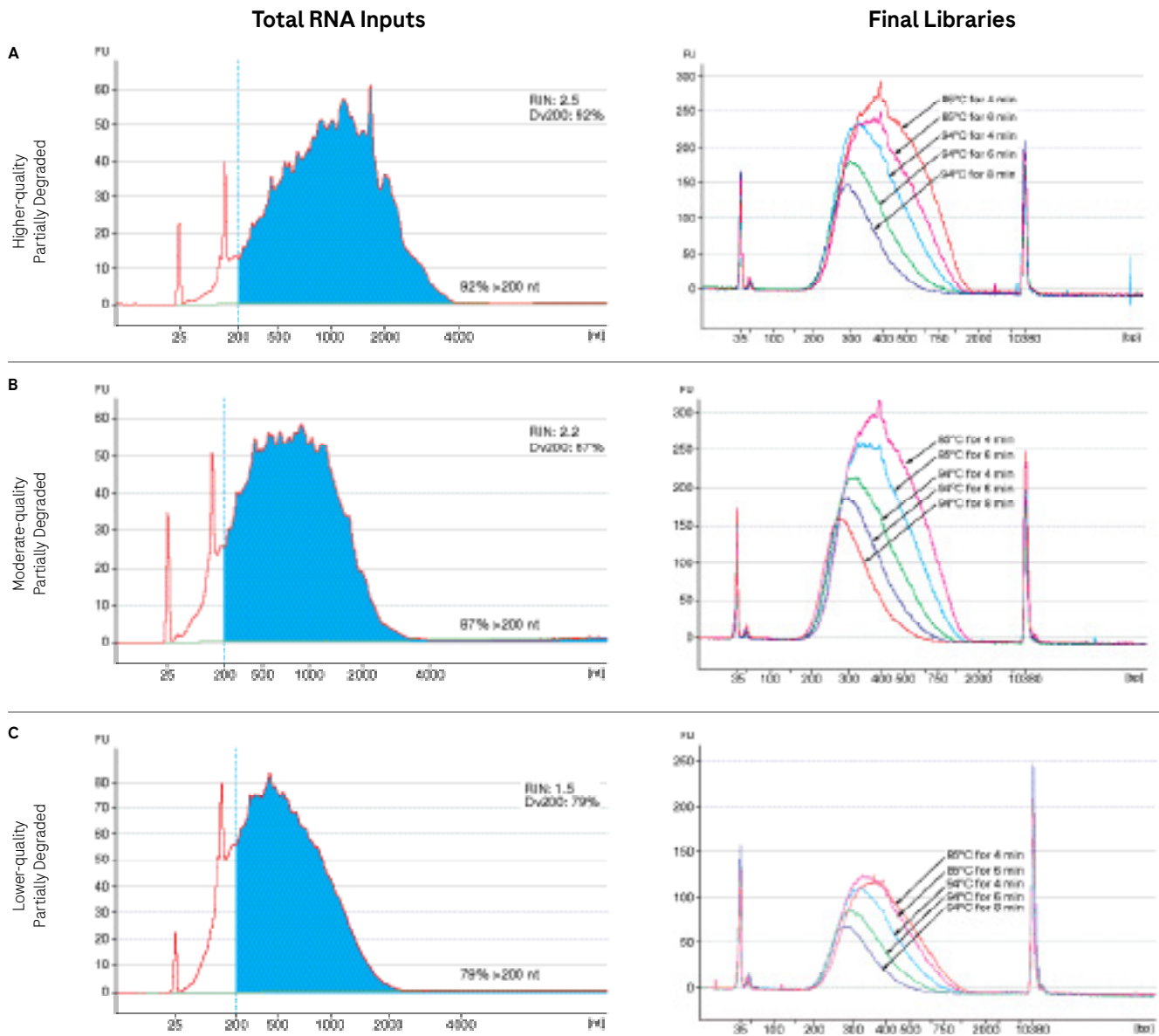

Fragmentation	Higher-quality Partially Degraded Final Library Size		Moderate-quality Partially Degraded Final Library Size		Lower-quality Partially Degraded Final Library Size	
	Mean (bp)	Mode (bp)	Mean (bp)	Mode (bp)	Mean (bp)	Mode (bp)
+94°C for 8 min	~320	~280	~310	~270	~310	~280
+94°C for 6 min	~350	~300	~340	~300	~330	~290
+94°C for 4 min	~380	~330	~380	~310	~360	~310
+85°C for 6 min	~420	~390	~410	~330	~390	~330
+85°C for 4 min	~450	~390	~450	~390	~410	~360

Figure 6. Input RNA and final library distributions for a range of partially-degraded sample qualities.

Libraries were constructed using 100 ng of chemically-degraded UHR RNA to target various qualities of partially degraded inputs, including higher-quality (A), moderate-quality (B), and lower-quality (C). As expected, the two commonly used RNA quality metrics, RIN and Dv_{200} decrease as RNA quality decreases. For all RNA qualities assessed, increased fragmentation time and/or temperature resulted in shorter, narrower distributions. Note that results may differ with other sample sources. Total RNA electropherograms were generated with an Agilent RNA 6000 Pico Kit, and final library electropherograms were generated with an Agilent 2100 High Sensitivity DNA Kit.

The effects of total RNA quality on final library distributions are illustrated in [Figure 7](#), and summarized below by fragmentation condition:

- +94°C for 8 min: Final library distributions were not overly impacted by input RNA quality, but a correlation between lower quality and reduced final library yield was apparent.
- +94°C for 4 min: While final library distribution modes were not overly impacted by input RNA quality, distributions became narrower, resulting in lower mean values, as RNA quality decreased.
- +85°C for 4 min: Both the mode and mean final library distribution metrics were impacted by RNA input quality, with both metrics decreasing as input quality decreased. Final library yield was not overly impacted until input quality dropped sufficiently low.

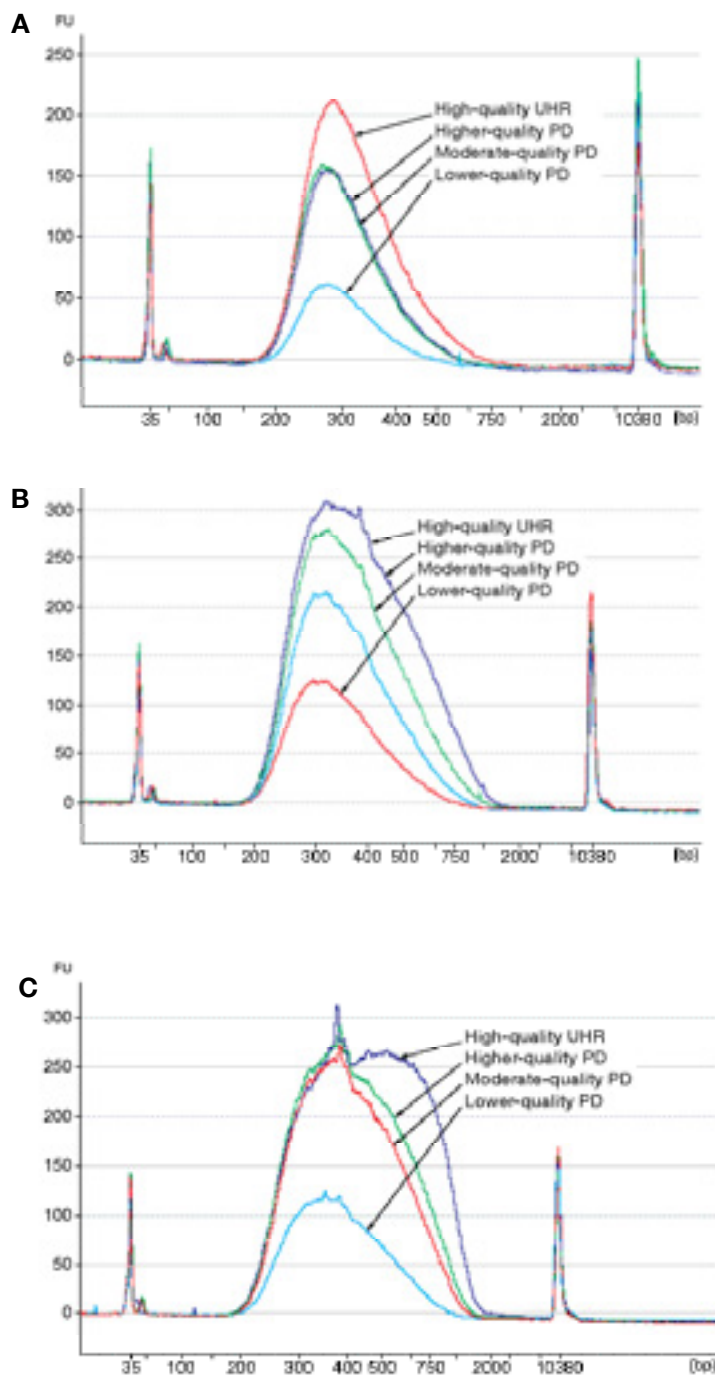


Figure 7. Assessment of the effects of total RNA quality on final library distributions.

Fragmentation conditions were selected to target a range of final library sizes: +94°C for 8 min (A), +94°C for 4 min (B), and +85°C for 4 min (C). Libraries were constructed using 100 ng of either intact UHR RNA (UHR) or partially-degraded (PD) UHR RNA. Results may differ with other sample sources. Electropherograms were generated with an Agilent High Sensitivity DNA Kit.



Appendix B. Limited Warranty

1. Limited Warranty

A. Products: Roche Sequencing Solutions, Inc. (“Roche”) warrants that its Products conform to its published specifications and are free from defects in material or workmanship. Customer’s sole and exclusive remedy (and Roche’s sole and exclusive liability) under this limited warranty shall be to either (a) replace the defective Products, or (b) provide Customer with a refund, as solely determined by Roche.

B. Under no circumstances shall Roche’s liability to Customer exceed the amount paid by Customer for the Services and Products to Roche. Roche will bear all reasonable shipping costs if service is re-performed at Roche or the Products are replaced. This warranty does not apply to any defect or nonconformance caused by (i) the failure by Customer to provide a suitable storage, use, or operating environment for the Materials or Customer’s submission of substandard quality Materials or contaminated or degraded Materials to Roche, (ii) Customer’s use of non-recommended reagents, (iii) Customer’s use of the Products, Materials or Data for a purpose or in a manner other than that for which they were designed, (iv) the failure by Customer to follow Roche’s published protocols; or (v) as a result of any other abuse, misuse or neglect of the Products, Materials or Data by Customer. This warranty applies only to Customer and not to third parties.

C. TO THE FULLEST EXTENT PERMITTED BY APPLICABLE LAW, ROCHE DISCLAIMS ALL OTHER REPRESENTATIONS, AND WARRANTIES, EXPRESS OR IMPLIED, WITH RESPECT TO THE PRODUCTS, SERVICES AND DATA, INCLUDING BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE OR NON-INFRINGEMENT. CUSTOMER’S SOLE REMEDY FOR BREACH OF WARRANTY IS STATED ABOVE.

D. Any action by Customer against Roche for Roche’s breach of this warranty must be commenced within 12 months following the date of such breach. Notwithstanding such 12-month period, within twenty (20) days of the delivery of Data and/or Products to Customer, Customer must notify Roche in writing of any nonconformity of the Services and Products, describing the nonconformity in detail; otherwise all Services and Products shall be conclusively deemed accepted without qualification.

2. Further Liability Limitation

TO THE FULLEST EXTENT PERMITTED UNDER APPLICABLE LAW, ROCHE SHALL NOT HAVE ANY LIABILITY FOR INCIDENTAL, COMPENSATORY, PUNITIVE, CONSEQUENTIAL, INDIRECT, SPECIAL OR OTHER SIMILAR DAMAGES, HOWEVER CAUSED AND REGARDLESS OF FORM OF ACTION WHETHER IN CONTRACT, TORT (INCLUDING NEGLIGENCE), STRICT PRODUCT LIABILITY OR OTHERWISE, EVEN IF ROCHE HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. CUSTOMER UNDERSTANDS THAT ANY RISKS OF LOSS HEREUNDER ARE REFLECTED IN THE PRICE OF THE SERVICES AND PRODUCTS AND THAT THESE TERMS WOULD HAVE BEEN DIFFERENT IF THERE HAD BEEN A DIFFERENT ALLOCATION OF RISK.

If you have any questions concerning service of this product, contact your local Roche Technical Support.

Go to sequencing.roche.com/support for contact information.

Evidence of original purchase is required. It is important to save your sales receipt or packaging slip to verify purchase.



KAPA and KAPA HYPERPLEX are trademarks of Roche. All other product names and trademarks are the property of their respective owners.

sequencing.roche.com

Roche Diagnostics Cape Town
19 Rubicon Boulevard, Brackengate 2, Brackenfell
7560, Cape Town, South Africa

For Research Use Only. Not for use in diagnostic procedures.